

Vehicle Cruise Control

Topics Covered

- Basic data gathering using **QUANSER RAPID CONTROL PROTOTYPING TOOLKIT®**
- PI speed control of a DC motor

Prerequisites

- The QUBE-Servo 2 has been setup and tested. See the QUBE-Servo 2 Quick Start Guide for details.
- Inertia disk load is on the QUBE-Servo 2.
- You are familiar with the basics of **LABVIEW™**.

Note: This workbook contains a single independent laboratory experiment as an introduction to the QDS platform. If you are interested in the complete QDS platform, please contact info@quanser.com

1 Background

1.1 Quanser® HIL Driving Simulator

The **QUANSER®** HIL Driving Simulator (QDS) is a modular and expandable **LABVIEW™** model of a car driving on a closed track. The model is intended as a platform for the development, implementation and evaluation of a variety of control systems. The QDS consists of a variety of components that are integrated together to create a representation of a vehicle being driven on a track. One possible configuration is shown in Figure 1.1. The model utilizes the **QUANSER RAPID CONTROL PROTOTYPING TOOLKIT®** (RCP) to facilitate hardware-in-the-loop interfacing (HIL). The **QUANSER®** 3D Viewer is also used to create an immersive visual environment for testing and evaluating controllers. Students are expected to observe and think critically about the effects of system parameters on not just the discrete plant, but the overall system.



Figure 1.1: **QUANSER®** Driving Simulator

Some examples of the real-world control problems that can be addressed using the QDS include parking assist systems, radar guided cruise control, active suspension, traction control and autonomous navigation.

1.2 Electronic Throttle Control

Electronic throttle control (ETC), traction control and cruise control have become standard features on modern cars. More recently, with the advent of radar-guided cruise control, and pre-crash systems electronic vehicle control units have begun to play an increasingly critical role in the real-time control of vehicle speed. Though the implementation of these systems can be complex, the essential system can be viewed as a closed-loop speed controller that sends throttle commands to the engine, shown in Figure 1.2.

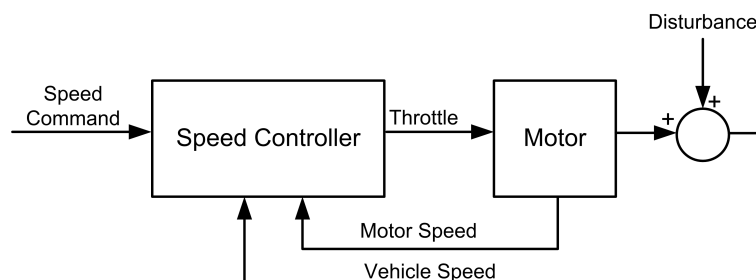


Figure 1.2: ETC controller

1.3 DC Motor Speed Control

For this laboratory, you will use the QDS in conjunction with a QUBE-Servo 2 to develop a Proportional-Integral (PI) speed controller to regulate the speed of the simulated vehicle. This exercise is analogous to developing an ETC for an electric vehicle, shown in Figure 1.3. In this case, the disturbance to the speed of the vehicle is the slope of the road, which is converted into a voltage and added to the control command. The controller that you will design takes the actual speed of the QUBE-Servo 2 motor as the derivative of the encoder measurements, and the commanded speed from the internal driver controller in the QDS. The PI controller then outputs the motor voltage V_m , which is added to the slope of the road and sent to the QUBE-Servo 2. The slope of the road is provided by the sensors inside the QDS.

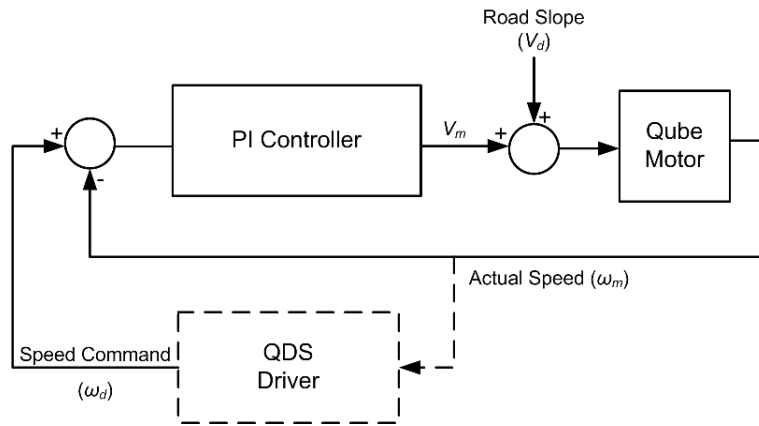


Figure 1.3: QDS speed controller

1.4 PI Control

The PI controller is one of the most common control algorithms. It combines the error reduction of proportional control, with the offset elimination of integral control. The proportional term tracks the instantaneous error, while the integral term controls the offset by tracking the total error over time. The longer there is error in the system, the more the integral term tries to compensate. Though the lack of a derivative term can cause overshoot and a longer settling time, for systems with simple dynamics the algorithm can provide near-optimal performance with no steady-state error.

In the time-domain, the linear behavior of a PI controller can be described by:

$$u(t) = k_p(r(t) - y(t)) + k_i \int_0^t (r(t) - y(t))dt \quad (1.1)$$

where $u(t)$ is the control signal, $r(t)$ is the reference, and $y(t)$ is the measured process output.

2 In-Lab Exercise

The Quanser Driving Simulator - Speed.vi VI shown in Figure 2.1 is used to implement the closed-loop PI speed controller on the QDS platform. The goal of this section is to observe the response of the actual system in the context of the driving simulator, and observe the effects of the PI controller on the performance of the vehicle.

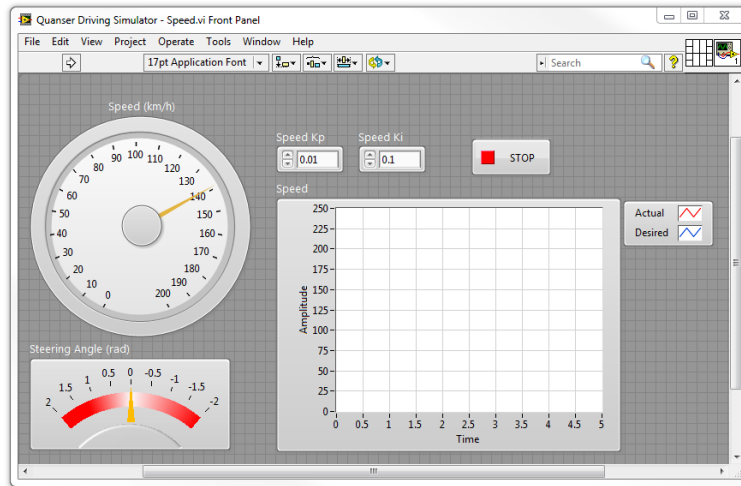


Figure 2.1: QUANSER® Driving Simulator with HIL Speed Control

Note: If you are using the QUBE-Servo 2 for NI myRIO, please begin by opening the Cruise Control - myRIO.lvproj LABVIEW™ project in the myRIO folder. The control subsystems are in the QUBE-Servo Speed Control - myRIO.vi VI deployed to the NI myRIO, and the simulation parameters in the Cruise Control.vi VI.

1. If the car has 18 in rims and 1 in tires, derive the value of *Car to Wheel* in the Speed Control subsystem block diagram. This gain converts the car speed in m/s to wheel rotational speed in rad/s.
2. Enter a k_p gain of 0.3 into the Kp gain control, and a value of 0 into the Ki gain control.
3. Run the VI.
4. Observe the performance of the car as it makes a lap of the track.
5. Set $k_p = 0.08$ and $k_i = 0.56$ and run the VI. What do you observe about the performance of the vehicle in both cases?
6. Is the controller able to track the desired speed effectively?
7. On the block diagram, change the magnitude of the Slope Gain constant to 0.
8. How well is the controller able to track the speed command if the disturbance is eliminated?
9. On the block diagram, replace the output of the Autodriver subsystem with a constant of 50 m/s.
10. Reset the Slope Gain block back to the original gain of 10.
11. How well is the controller able to compensate for the slope disturbance? How does this compare to the setpoint tracking performance?
12. If necessary, re-tune the controller gains to achieve the desired performance.
13. Record the final control gains and response plots.
14. Click on the Stop button to stop the VI.
15. Power OFF the QUBE-Servo 2.

© 2016 Quanser Inc., All rights reserved.

Quanser Inc.
119 Spy Court
Markham, Ontario
L3R 5H6
Canada
info@quanser.com
Phone: 1-905-940-3575
Fax: 1-905-940-3576

Printed in Markham, Ontario.

For more information on the solutions Quanser Inc. offers, please visit the web site at:
<http://www.quanser.com>

This document and the software described in it are provided subject to a license agreement. Neither the software nor this document may be used or copied except as specified under the terms of that license agreement. Quanser Inc. grants the following rights: a) The right to reproduce the work, to incorporate the work into one or more collections, and to reproduce the work as incorporated in the collections, b) to create and reproduce adaptations provided reasonable steps are taken to clearly identify the changes that were made to the original work, c) to distribute and publically perform the work including as incorporated in collections, and d) to distribute and publicly perform adaptations. The above rights may be exercised in all media and formats whether now known or hereafter devised. These rights are granted subject to and limited by the following restrictions: a) You may not exercise any of the rights granted to You in above in any manner that is primarily intended for or directed toward commercial advantage or private monetary compensation, and b) You must keep intact all copyright notices for the Work and provide the name Quanser Inc. for attribution. These restrictions may not be waved without express prior written permission of Quanser Inc.